


THE 3001KL-AM

ANTIMICROBIAL DIGITAL PLATFORM SCALE

The all-digital 3001KL-AM supports up to 1000 pounds and provides the optimum quality, reliability, and ease of use with a state-of-the-art digital height rod. Designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities, The scale features a large non-skid low profile platform, "live" handrails for patient stability, and rugged scuff-proof wheels that enable portability. The scale also offers an antimicrobial keypad, handrails, and transport handle to give patients and providers additional safety. The scale's high-contrast color screen provides a user-friendly interface that includes BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently

Antimicrobial

Keypad
Handrails
Transport Handle

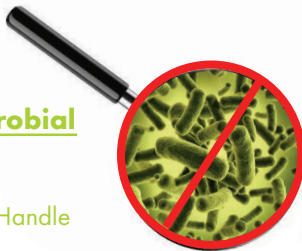

The 3001KL-AM offers two optional height rod options with quick and easy to assemble brackets, a digital height rod 245EHR-3001 or a mechanical height rod 3001HR. The 245EHR-3001 market-leading digital height rod provides the most accurate height measurement in its class, with a graduation/resolution of 1mm/0.1" and a full measuring range of 11 3/4" – 79 1/2" / 30cm – 202cm and a digital range of 3ft 11.2" – 6ft 7.5" / 120 cm - 202 cm. The 3001HR features a durable height rod construction and an easy-slide head-piece that allows caregivers to quickly and effortlessly measure and read a patient's height. The 3001HR provides a measuring range of 2"-90" in 1/16" increments (6 cm-230cm in 1 mm increments).

The scale provides easy integration and seamless communication with today's EMR applications. The 3001KL-AM includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 3001KL-AM is your choice for safety and reliability.

120+ LOCAL
FIELD SPECIALISTS

100%
CONNECTIVITY

#1 TRUSTED
BRAND

100 YEARS
OF SERVICE



Health o meter®

Professional

Health o meter® Professional 3001KL-AM

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

1000 lb / 454 kg

Resolution

0.2 lb / 0.1 kg

EMR Connectivity

Via USB, optional Pelstar® wireless technology

Display

Digit Size: 3/8" / 10 mm
Screen Size: 2 3/4" x 1 1/2" High-Contrast
Color TFT-LCD Screen

Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

Functions

LB/KG conversion
LB/KG lock out
BMI
Zero
Tare
Pre-tare
Hold/release
Reweigh
Recall
Auto zero
Auto off
Time/date
Variable auto off time
Audible/mute sound option
Parent/child weighing function
Everlock® - permanent weight lock

Platform Size (w x d x h)

22" x 15 3/4" x 2 5/8"
(559 mm x 400 mm x 64 mm)

Product Footprint (w x d x h)

22" x 33 1/8" x 50"
(559 mm x 841 mm x 1270 mm)

Product Weight

46 lb / 21 kg

Portability

2 wheels for ease of movement

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

MODEL	DESCRIPTION
3001KL-AM	Antimicrobial digital platform scale, LB & KG
3001KG-AM	Antimicrobial digital platform scale, KG only
3001KL-AM-BT	3001KL-AM with built-in Pelstar® wireless technology
3001KG-AM-BT	3001KG-AM with built-in Pelstar® wireless technology
3001KL-AMHR	3001KL-AM with mechanical height rod
3001KG-AMHR	3001KG-AM with mechanical height rod
3001KL-AMHR-BT	3001KL-AMHR with built-in Pelstar wireless technology
245EHR-3001	Optional digital height rod for 3001 models*
3001HR	Optional mechanical height rod for 3001 models
SS-3001KL-AM	ScaleSurance 2 year extended warranty (3001 models)**
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit homscales.com/ConnectivitySolutions

Optional Height Rods Available



To learn more, call **1-800-815-6615**
or visit homscales.com/3001KL-AM

Health o meter®
Professional



*Compatible with 3001 series scales including a date code with the last two digits of 21 or higher. An upgrade option is available for models including a date code with the last two digits of 17 or higher. Not available for AMX series scales. For assistance, contact Customer Service at 800-815-6615.

** Not available in all countries