

THE 500KLHB

PRIMARY CARE DIGITAL SCALE WITH LIVE HANDLEBARS

The number-one-selling physician scale in North America, the all-digital Health o meter® 500KLHB primary care scale provides the optimum in quality, reliability, and ease of use. The unique live handlebars provide stability that many patients need, without affecting the accuracy of the weight reading. Designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities, the scale features a built-in height rod in the scale column, a large non-skid platform, and reinforced base for added stability. The scale gives providers and patients an eye-level weight indicator for patients up to 660 pounds, with a user-friendly digital interface that includes BMI, hold/release, zero, tare, reweigh, recall, KG/LB, clear, and enter functions. The 500KL offers the best in its class extra-wide 8" wheelbase making transport easy. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

The first healthcare scale to provide easy integration and seamless communication with today's EMR applications, the 500KLHB includes the option to add wireless connectivity to interface with a computer, EMR software, or other electronic device.

With its low platform height, the 500KLHB is easy to set up and offers patients a more comfortable and stable weighing experience.

THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 500KLHB is the most widely used scale in physician offices nationwide.

660 LB
CAPACITY


150,000+
INSTALLED
NATIONWIDE

120+ LOCAL
FIELD SPECIALISTS

8" WIDE WHEELBASE
BEST IN ITS CLASS

#1 TRUSTED
PHYSICIAN SCALE

100 YEARS
OF SERVICE

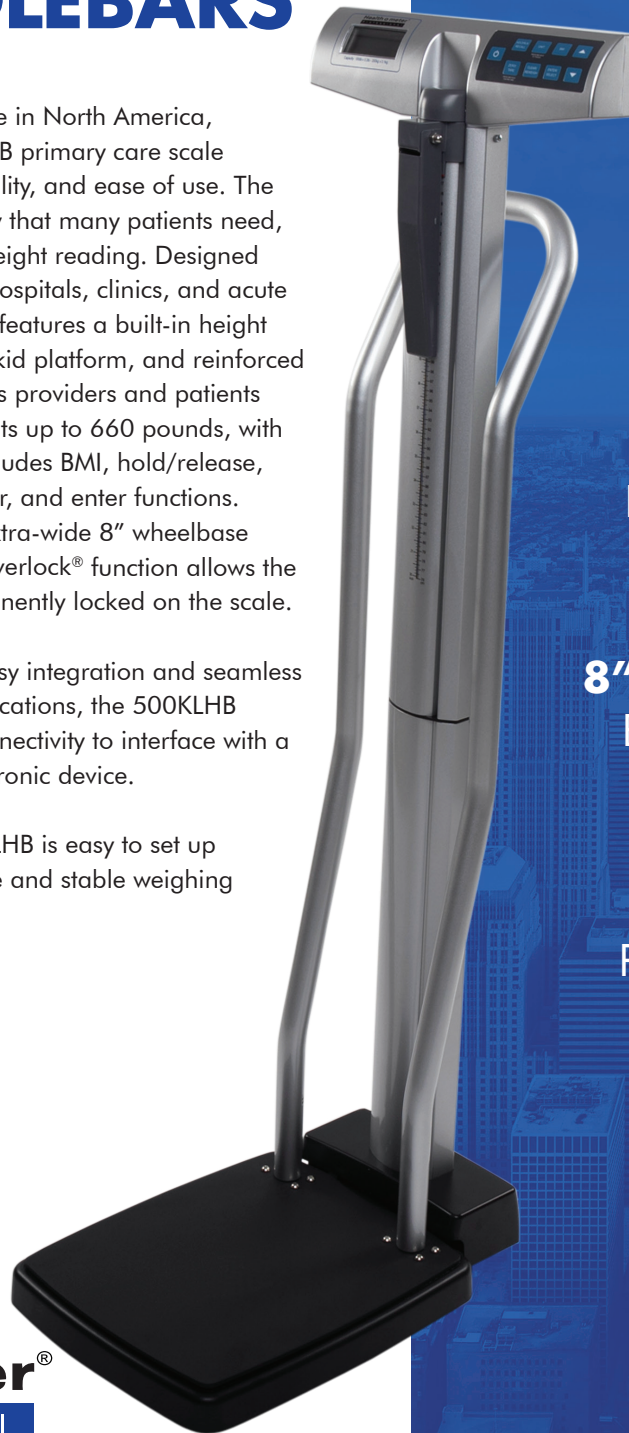

Health o meter®
Professional

Health o meter® Professional 500KLHB

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

660 lb / 300 kg

Resolution

0.2 lb / 0.1 kg

EMR connectivity

Via optional Pelstar® wireless technology

Display

1" / 25 mm, LCD display

Batteries

6 AA batteries included, optional adapter (order# ADPT31)

Functions

LB/KG conversion
LB/KG lock out
BMI
Zero
Tare, Pre-tare
Hold/Release
Auto Zero
Auto Off
Reweigh
Recall
Parent/Child weighing function
Everlock® - permanent unit lock

Platform Size (w x d x h)

13 3/4" x 16 1/2" x 2 3/8"
(349 mm x 419 mm x 60 mm)
Dent and rust proof platform

Product Footprint

14 3/8" x 21 1/8" x 53 5/8"
(359 mm x 537 mm x 1362 mm)

Product Weight

31 lb / 15 kg

Height Rod

30"-84 1/4" / 76 cm-214 cm
Graduation: 1/8" / 1 mm

Portability

2 wheels, wide 8" wheelbase for ease of movement

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

MODEL	DESCRIPTION
500KLHB	500KL with handlebar accessory
SS-500KL	ScaleSurance® 2 year extended warranty for 500 series*
ADPT31	Optional adapter for 500KLHB
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit hom scales.com/ConnectivitySolutions

The handlebar is easy to attach to any new scale using the pre-drilled mounting holes and can also be attached to any existing scale using the easy-to-use drill guides.



Health o meter®
Professional



*Not available in all countries.

To learn more, call **1-800-815-6615**
or visit hom scales.com/500KLHB